

# Heyer Auto CPAP Devices



**TADAWI**  
LIFE TECH

Your Life...Our Commitment

كياتكم..التزامنا

-  Muscat, Oman
-  +968-91235160, 79012020
-  info@tadawilifetech.com
-  Tadawilifetech.com

# Product Specification

Product Model	20C	20A	20A Plus	20A Pro
Therapy Mode	CPAP	Auto CPAP/CPAP	Auto CPAP/CPAP	Auto CPAP/CPAP
Pressure Range (cmH2O)	4~20	4~20	4~20	4~20
Auto Ramp Time (min.)	0~60	/	0~60	0~60
Auto As-Elx Pressure Release	1~7Level	1~7Level	1~7Level	1~7Level
Identification of Central Apnea	Yes	Yes	Yes	Yes
ASDry Technology	Yes	Yes	Yes	Yes
Auto On/Off	Yes	Yes	Yes	Yes
Heating Tubing	Yes(Optional)			
Cellular	LTE Cat M1 & Cat NB2			
Bluetooth	Yes(Optional)			
Leaking Compensation	Yes	Yes	Yes	Yes
Titration	/	/	/	Yes
Therapy Data Storage	Yes, SDCard			
Automatic Humidification		/	Yes	Yes
Screen Size	Yes	3.5 Inches LCD Screen		
Machine Size	267 mm * 157 mm * 106 mm			
Weight	1.5kg			
Working Noise	28±2 dB			
Power Adapter	AC Input: 100-240V, DC Output: 24V			



# CPAP Device

—Continuous Positive Airway Pressure Machine  
Accurate Intelligent Control for Every Comfortable Breath



**Wide-viewing Operation Interface**  
3.5 inches LCD Color Screen



**One-click Initiation and Cessation**  
Click Standby/Home button to Start/Stop Therapy

**Updated Water Tank**  
Big Opening and Easy for Refill and Cleaning  
Preheating Time up to 30 Min.




**Noise Level**  
28dB

**Efficient Communication System**  
Cellular/Bluetooth Module/SD Card




**Multiple Choices of Tubing**  
5 Heating Levels  
Heating Tube with Dia. at 15mm  
Standard Tube with Dia. at 15mm or 19mm


- Adapt to a variety of masks
- Heating Tube Available

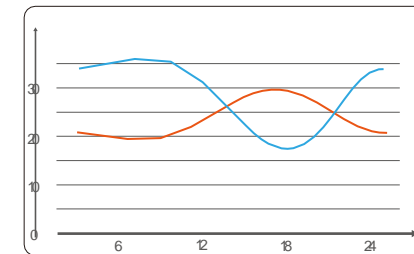

**Intelligent Temperature and Humidity Control**  
The CPAP system intelligently adjusts temperature and humidity levels in response to environmental variations, ensuring optimal comfort during breathing.

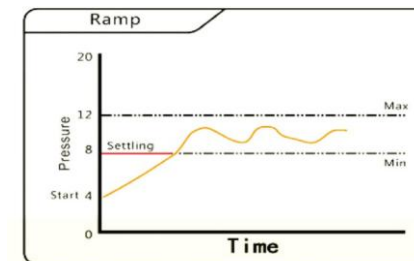

**Auto Ramp Time: Max. Up to 60min**  
The system intelligently detects the user's sleeping status and automatically adjusts the ramp time to align with their duration in bed, ensuring precise treatment pressure once asleep.

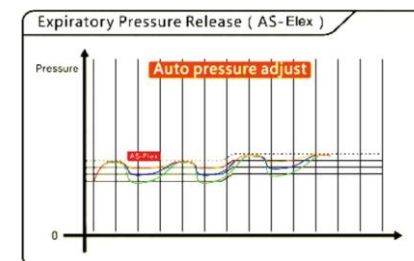

**As-Elex Pressure Release: 7 Grades x 0.5 cmH2O :**  
AS-Elex provides 7 grades of pressure release during expiration, each grade with 0.5cmH2O differences. Users have more choice to adjust the pressure release to avoid discomfort, which effectively increases patient compliance.